

# Transcatheter interventions for tricuspid regurgitation - heterotopic technology: TricValve

Hans R. Figulla<sup>1\*</sup>, MD; Katharina Kiss<sup>2</sup>, MD; Alexander Lauten<sup>3</sup>, MD

1. Kerckhoff Heart Center; Bad Nauheim, Germany; 2. Intern Kardiologische Angiologische Gruppenpraxis, Vienna, Austria;  
3. Department of Cardiology, Charité-Universitätsmedizin Berlin, Campus Benjamin Franklin, Berlin, Germany

This paper also includes supplementary data published online at: [http://www.pcronline.com/eurointervention/Y\\_issue/??](http://www.pcronline.com/eurointervention/Y_issue/??)

## Introduction

Considering the incidence of tricuspid regurgitation (TR) and its implications on functional status and long-term survival, TR is currently undertreated. With the increasing adoption of catheter-based treatments of other types of structural heart disease there is a growing interest and need for effective interventional treatment also for tricuspid regurgitation. Therefore, multiple treatment approaches are under investigation; however, in this chapter we focus exclusively upon tricuspid replacement.

## Transcatheter replacement of the tricuspid valve

Due to the high operative mortality of isolated tricuspid surgery, it is probably correct to predict that the percutaneous approach will be the treatment of choice in patients requiring tricuspid surgery. Considering the large unmet need for an interventional option for inoperable patients with tricuspid regurgitation, percutaneous valve implantation could be a feasible alternative, besides interventional valve repair. From the interventional

perspective, there are two basic approaches, depending on the anatomic site of prosthetic valve implantation – an orthotopic versus a heterotopic valve replacement.

In orthotopic valve replacement, the prosthetic valve is deployed at the level of the TV annulus between the right ventricle (RV) and right atrium (RA). This approach was investigated by Boudjemline et al by means of implanting a double-disc nitinol stent with a semi-lunar valve into the tricuspid annulus. Due to the anatomic structure and the flexibility of the surrounding myocardium, this site of implantation offers little resistance for orthotopic long-term fixation of stent-based valves. In functional TR, annulus diameter may reach >50 mm associated with a loss of anatomic landmarks between the RV and RA. A device intended for orthotopic TV replacement would therefore require unique solutions for stent-valve fixation as well as tissue valve engineering (e.g., a 50 mm tissue valve would require a profile height of >30 mm). To our knowledge only one company is presently working on such a demanding system for orthotopic valve replacement (TRicares GmbH, Aschheim, Germany).

\*Corresponding author: In den Bornwiesen 4, D-07749 Jena, Germany.

E-mail: [figulla@figulla.org](mailto:figulla@figulla.org)

54 So far in patients, orthotopic tricuspid valve replacement has  
 55 been used with promising results only as a valve-in-valve or  
 56 valve-in-ring procedure using either the balloon-expandable  
 57 Edwards SAPIEN (Edwards Lifesciences, Irvine, CA, USA) or  
 58 the Melody® (Medtronic, Minneapolis, MN, USA) valve. As with  
 59 most low-volume procedures, there are no published data from  
 60 larger patient series regarding the outcome and long-term function  
 61 of these devices in the tricuspid position.

62 A feasible alternative to orthotopic approaches is caval valve  
 63 implantation (CAVI), which involves the implantation of stent  
 64 valves in a heterotopic position into the inferior vena cava (IVC)  
 65 and superior vena cava (SVC). This concept has been investi-  
 66 gated preclinically with encouraging results and has been used  
 67 on compassionate grounds in humans. However, pulsatile sys-  
 68 tolic blood flow reversal in the caval veins is a prerequisite  
 69 for the proper function of the caval valves. This is reflected by  
 70 a profound V-wave atrial flow pattern which is required prior to  
 71 heterotopic implantation.

### 72 Heterotopic valve design

73 The CAVI procedure is technically simple and can be performed  
 74 rapidly without interference to cardiac structures. Devices do not  
 75 interfere with any pre-existing transtricuspid pacemaker or defi-  
 76 brillator leads, which might represent a limitation for orthotopic  
 77 procedures on the tricuspid valve. Under the condition of severe  
 78 chronic TR, there is a considerable variation of the anatomic

diameter of the SVC and IVC, which frequently exceeds the suit-  
 able range for implantation of current, commercially available  
 devices. The SVC might dilate up to 40 mm in diameter and the  
 IVC up to 35 mm. With the Edwards SAPIEN XT valve, these  
 challenges have been partially compensated for by pre-stenting  
 the caval veins prior to valve deployment for downsizing and  
 improved valve anchoring by Laule et al who showed haemo-  
 dynamic improvement even with a single valve only in the IVC.


In the caval position, self-expandable devices specifically  
 designed for caval valve implantation with little radial force are  
 likely to be a superior alternative as they do not require pre-  
 stenting of the landing zone. Self-expandable valves for commer-  
 cial use are currently under development by the company P & F  
 Products & Features Vertriebs GmbH, Vienna, Austria. The ana-  
 tomic variation is addressed for the SVC by a device with a one-  
 size valve diameter of 30 mm in the tubular part of the stent, but  
 with a variable hip protrusion of up to 45 mm. The IVC device  
 is deployed with the jacketed valve part in the right atrium and  
 a waist in the hiatus of the diaphragm, allowing hepatic veins to  
 drain unobstructed into the IVC (**Figure 1**).

### 79 Haemodynamics and anatomic consequences 80 following bicaval valve implantation (CAVI)

81 Following implantation of the TricValves, the valves resemble car-  
 82 diac valves (**Moving image 1**). Caval V-waves are flattened, while  
 83 atrial V-waves increase following the reduced capacity volume for



84  
85  
86  
87  
88  
89  
90  
91  
92  
93  
94  
95  
96  
97  
98  
99  
100  
101  
102  
103  
104  
105  
106 **Figure 1.** The TricValve (produced by P & F Products & Features Vertriebs GmbH, Weßling, Germany), under development.

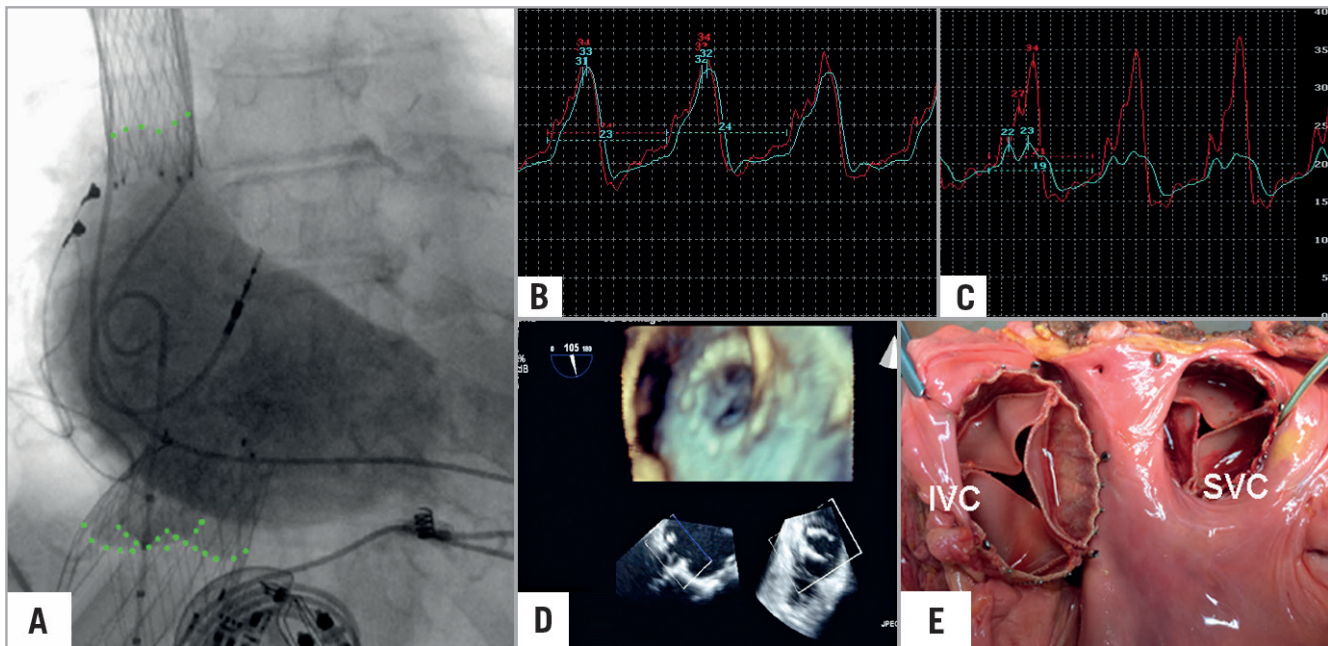

**Figure 2.** Caval valve implantation. *A)* Position of self-expanding valves in the SVC and the IVC. Note the valve leaflets marked by angiography. *B)* & *C)* Pressure measurement confirms a reduction of V-wave and mean pressure in the inferior vena cava from 32 mmHg to 23 mmHg and 24 mmHg to 19 mmHg, respectively. *D)* Device function is observed in transoesophageal echocardiography. *E)* Macroscopic specimen demonstrating the device position in the SVC and IVC. (image copyright European Heart Journal)

the regurgitant blood flow. In the acute phase, mean caval pressures are only slightly reduced (2-3 mmHg). In the following weeks, the remodelling of the RV due to the reduced volume load initiates a decrease of RA and caval vein pressures. Post-mortem examinations showed an ingrowth of the stent frame, while a thrombus on the valve was never seen. However, all patients were on oral anticoagulation.

The persistence of right atrial volume overload, the ventricularisation of the RA and the increase of RV afterload are potential limitations of the procedure. Its long-term impact on RA and RV function is currently unknown. Post-implant device function has been observed up to 24 months following implantation in one patient, resulting in a profound clinical improvement, from NYHA Class IV to II and a normalisation of impaired liver function.

## Conclusion

CAVI with the TricValve is a relatively simple procedure. However, valve design has to cover a great range of caval vein anatomy. The haemodynamic concept is convincing and allows the RV to recover. Clinical experience is presently restricted to compassionate cases; commercialisation is on its way.

## Conflict of interest statement

H.R. Figulla, K. Kiss, and A. Lauten receive consultancy fees from P&F Products + Features, Vienna, Austria.

## Supplementary data

**Moving image 1.** An echo view of the IVC with a TricValve. Note the movement of the leaflets, resembling a cardiac valve. Colour code demonstrates no regurgitant flow.